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NEOGEN

# NeoSal®

## *Simplifying Oral Fluid Collection*

NEOGEN® is proud to announce the release of our NeoSal® Oral Fluid Collection System. Our oral fluid collection system is designed for maximum ease of use at the collection site and testing laboratory.

### Confident Collection

- Sample volume adequacy indicator: Ensures appropriate and consistent sample volume is collected.
- Stimulant-free collection pad: Collects sufficient volume without artificial stimulants.
- Quick collection: Typical collection time is 1–2 minutes.

### Engineered Ease of Use

- Unique collector handle design: Handle remains outside of the buffer tube, reducing risk of sample contamination and forming a leak-proof seal for long term storage.
- No additional processing: No centrifugation or plunging necessary, discard collection pad and process sample.
- Automation friendly: Buffer tube fits onto most standard automated instrument sample racks (12 x 75), eliminating transfer to additional tube.
- Additional cap: Designed for long-term storage, stores in bottom of tube when not in use.

### Optimized Performance

- Oral fluid ELISA kits: Optimized with the NeoSal, samples can be assayed without additional sample preparation.
- Drug recovery: Innovative proprietary buffer extracts adequate analyte concentration for consistent detection.
- Sample stability: Buffer contains microbial growth inhibitors and ensures stability of the sample during transport.
- Compatibility: NeoSal buffer is compatible with liquid chromatographic separation approaches and mass spectrometer detection methods.

### Oral Fluid Collection System Ordering Information

| Product | Available Sizes | Product No. |
|---|---|---|
| NeoSal Oral Fluid Collection System | 25 pack | 128101-25 |
| NeoSal Buffer | 500 mL<br>1 L | 120176<br>120177 |
| Buffer Tube caps | 25 pack | 128102-25 |

**Call 800.477.8201 to order or visit NEOGEN.com**


©NEOGEN Corporation, 2021. NEOGEN and NeoSal are registered trademarks of NEOGEN Corporation, Lansing, MI 48912 U.S.


*Simplifying Oral Fluid Collection*

## Simple Collection Process


The NeoSal Oral Fluid Collection System is designed for maximum ease of use at the collection site and testing laboratory. Collection is completed with these steps:

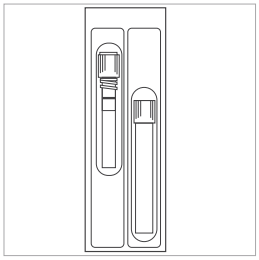
Peel open package and remove collector. Do not eat or drink at least 15 minutes prior to collection.

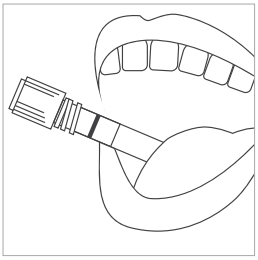
Pool saliva in mouth. Insert collector pad in mouth between cheek and gum.

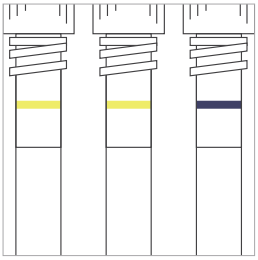When a distinct blue line forms, remove collector from mouth. Typical collection time is 1–2 minutes.

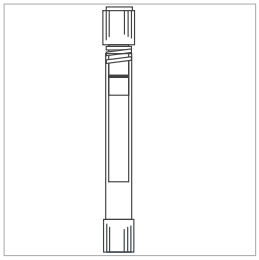Remove blue cap and snap into bottom of the tube. Insert saturated collector pad into tube.

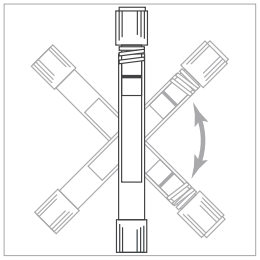Tighten white cap. Gently invert tube 4 or 5 times. Send sample to testing laboratory.

## NEOGEN's Drug Detection (Ready to Use) Forensic Assays

NEOGEN's forensic drug detection immunoassays include more than 100 ELISA test kits to screen over 300 drugs and metabolites in a wide range of forensic samples including whole blood, urine, oral fluid, and others.

| ELISA Test Kits | Product No. |
|---|---|
| 6-AM | 134019 |
| Acetaminophen | 132419 |
| Amphetamine Specific-2 | 132319-2 |
| Amphetamine Ultra | 130819 |
| Barbiturate Group | 130619 |
| Benzodiazepine Group (oxazepam/clonazepam) | 130119 |
| Buprenorphine | 131919 |
| Carisoprodol | 132519 |
| Citalopram | 132719 |
| Cocaine/Benzoylecgonine-2 | 130319-2 |
| Dextromethorphan | 131419 |
| Fentanyl | 131519 |
| Hydromorphone | 132219 |
| Ketamine | 131719 |
| LSD | 130219 |
| Meperidine | 131219 |
| Methadone | 131619 |
| Methadone/LAAM | 132919 |
| Methamphetamine/MDMA | 130919 |
| Methylphenidate/Ritalinic Acid | 134219 |
| Naltrexone/Nalbuphine | 133019 |
| Opiate Group | 130419 |
| Oxycodone/Oxymorphone | 130719 |

| ELISA Test Kits | Product No. |
|---|---|
| Phencyclidine (PCP) | 130519 |
| Phenytoin | 132119 |
| Propoxyphene | 131119 |
| Salicylates | 133619 |
| Sertraline | 131319 |
| Synthetic Cannabinoids (JWH-018) | 133519 |
| Synthetic Cannabinoids (JWH-250) | 133719 |
| Synthetic Cannabinoids (UR-144/XLR-11) | 133919 |
| Synthetic Cathinones (Methcathinone) | 181819 |
| THC | 131019 |
| Tramadol | 131819 |
| Trazodone | 132819 |
| Tricyclics Group | 132019 |
| Zolpidem | 132619 |

| Oral Fluid ELISA Kits | Product No. |
|---|---|
| Cocaine/BZE | 120119 |
| Methamphetamine/MDMA | 120219 |
| Opiates | 120319 |
| PCP | 120419 |
| THC | 120519 |
| Oral Fluid Multi-Analyte Calibrator Set | 128100 |

See additional kits at NEOGEN.com

**FORENSIC USE ONLY**

**Call 800.477.8201 to order or visit NEOGEN.com**

NEOGEN®
NEOGEN.com

©NEOGEN Corporation, 2021. NEOGEN and NeoSal are registered trademarks of NEOGEN Corporation, Lansing, MI 48912 U.S.