SNAPSHOOT
ENGINEERED FOR A SAFER TOMORROW

현장을 이해하고, 필요한 솔루션을 연결합니다.

산업안전 · 추락방지 시스템
산업용 소재 · 건축자재
제품 제안부터 공급·시공·설치·교육까지



COMPANY INFORMATION

스냅슛 (SNAPSHOOT)

ADDRESS (15431) 경기도 안산시 단원구 풍전로 37-9 트리플렉스 301동 118호
TEL 0507-1374-9229
FAX 0504-481-6749
E-MAIL bkkim@snapshoot.co.kr

TECHNICAL DOCUMENTS

인증서 · 데이터시트 · 기술자료

제품의 규격과 성능, 인증 및 적용 조건을 확인하기 위한 자료입니다.

ABOUT SNAPSHOOT

제품 공급을 넘어, 현장에 맞는 해결책으로

스냅슛은 산업현장의 작업 환경과 요구 조건을 바탕으로
적합한 제품과 시스템을 제안합니다. 현장 실사를 통해 설치 조건을 확인하고,
공급부터 시공·설치 및 사용 교육까지 필요한 과정을 함께합니다.

추락방지 시스템

안전대부착설비 · 수평·수직 라이프라인
레일 시스템 · 리지드 포스트

밀폐공간 진입·구조 솔루션

데빗암 · 삼각대 · 구조용 윈치
윈치 안전블록 · 구조장비

산업안전 보호구

전신형 안전대 · 랜야드 · 안전블록
작업 환경에 맞춘 개인보호구

산업용 소재·건축자재

산업용 테이프 · 접착제 · 기능성 소재
구조용 스크류 · 목조건축 자재

SNAPSHOOT
ENGINEERED FOR A SAFER TOMORROW

홈페이지 바로가기 →
문서 라이브러리 바로가기 →
시공사례 바로가기 →
문의하기 →

본 페이지는 스냅슛의 자료 제공 안내용 표지입니다. 첨부된 인증서 및 기술자료의 발행 주체, 인증 대상과 적용 범위는 해당 원문을 기준으로 확인해 주세요.



NEOGEN®

ANSR®

Software Upgrade Quick Reference

NEOGEN®

ANSR Desktop Software / Instrument Application Update

Summary

A software update for ANSR is available. This will upgrade ANSR systems to v2.1.0.4.

Software enhancements include improved sync functionality and system performance improvements.


The installer for the update can be downloaded at the FTP below:

ftp://ansrguest:Guest@ftp.neogen.com/ANSRUpgrade.zip.


A USB jump drive with installer may also be secured from Neogen. Reference SKU ANSR USB.

Desktop Software Installation


1. From the FTP or USB drive, copy the folder **ANSR v2.1.04** to the local desktop. Extract the files and save to local desktop.
2. Double click to run the file **T16-Desktop-installer_v2.1.0.4-447**.

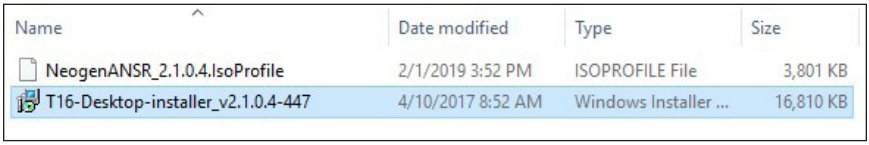


3. Select **Next**.

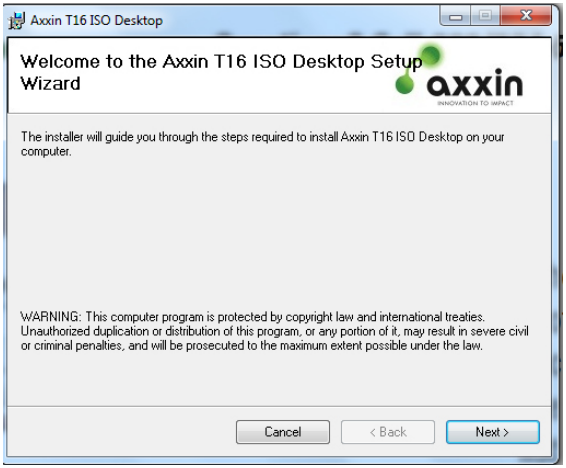


4. Specify the installation location. Then select **Next** to begin the installation.

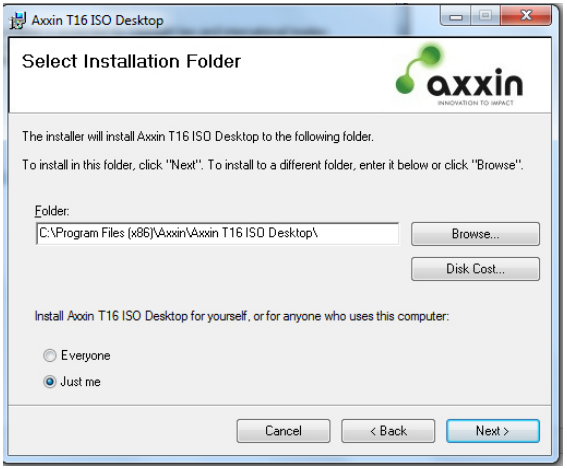

5. Select **Next** again to confirm the installation.

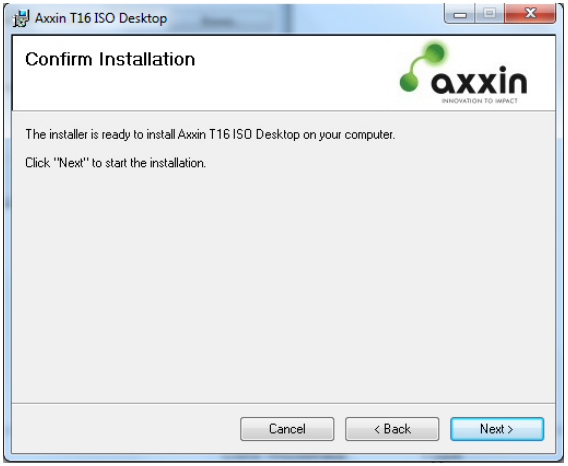


6. Read and accept the terms and conditions by clicking **I Agree**, and then **Next**.

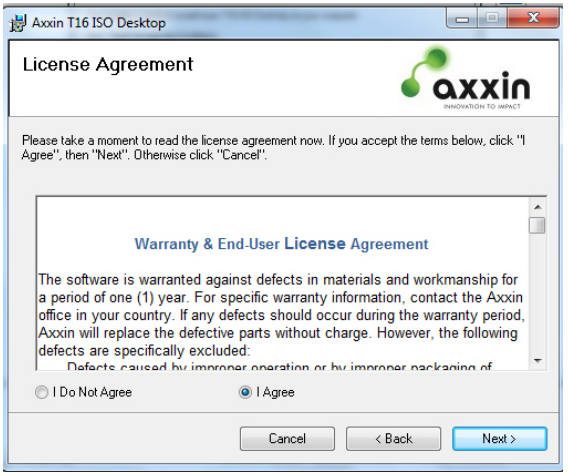


7. Installation will begin. Click **Next** once prompted.





8. To specify the customization file, click **Browse** and navigate to the desktop folder.

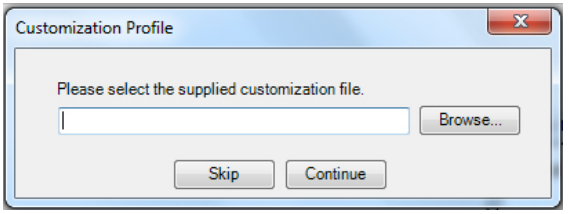

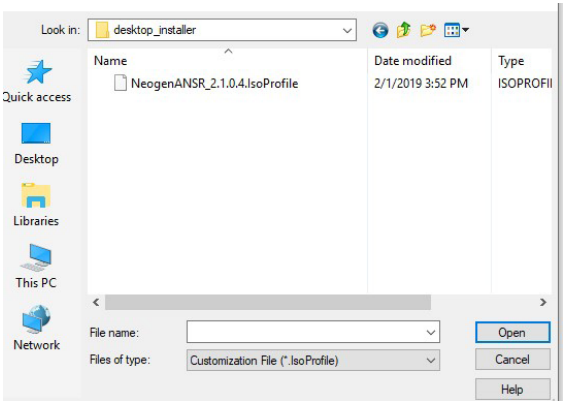



9. Follow the prompts on-screen to install the bridge driver.

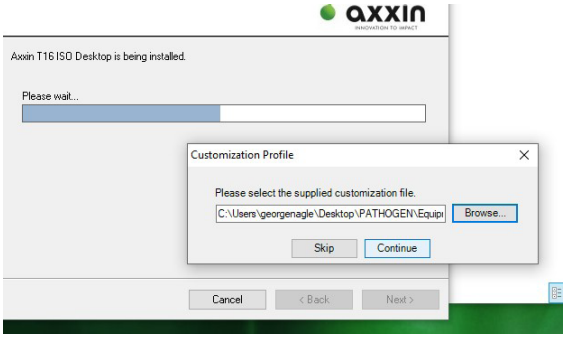


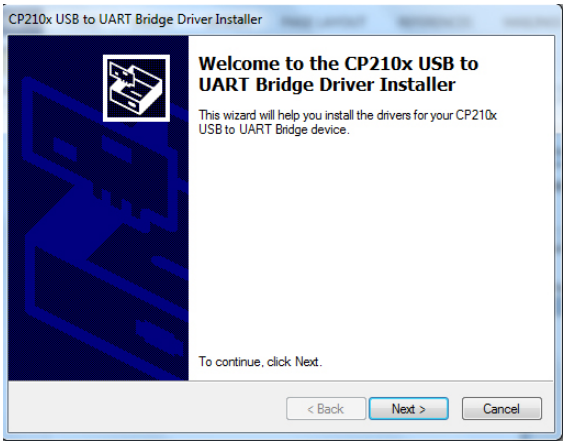

ANSR Desktop Software / Instrument Application Update

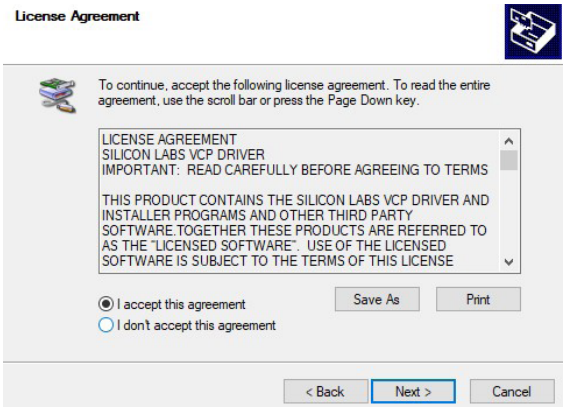


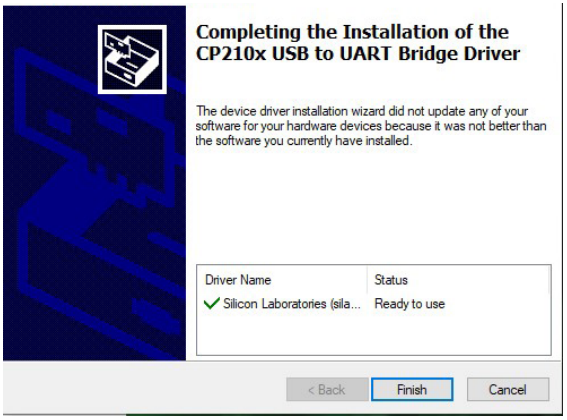

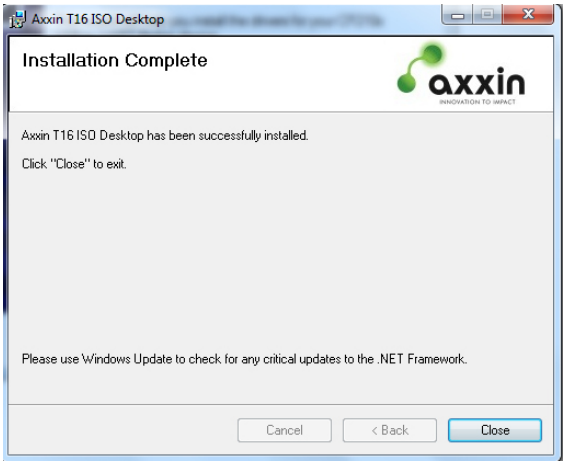


10. Locate the file **Iso.Desktop.System.T16** from **C:\Program Files (x86)\Axxin\Axxin T16 ISO Desktop**. For future use, create a shortcut to the desktop. Double click to open the desktop application.

11. Once the desktop application is open, select **Support** from the main menu.

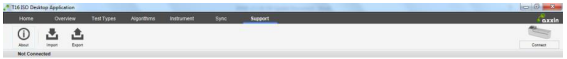

Refer to user manuals and the quick guides for this product
or contact support@axxin.com.

12. Select **Import** and specify the Open the .tar file **ANSR 20180315 .tar**. This completes the desktop application upgrade.

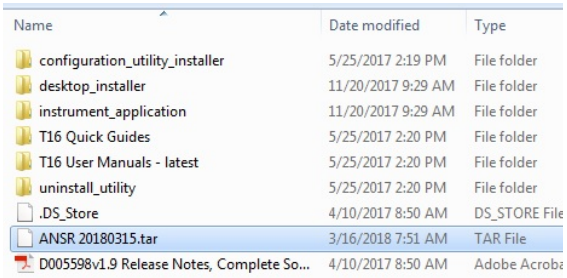

Name	Date modified	Type
configuration_utility_installer	5/25/2017 2:19 PM	File folder
desktop_installer	11/20/2017 9:29 AM	File folder
instrument_application	11/20/2017 9:29 AM	File folder
T16 Quick Guides	5/25/2017 2:20 PM	File folder
T16 User Manuals - latest	5/25/2017 2:20 PM	File folder
uninstall_utility	5/25/2017 2:20 PM	File folder
.DS_Store	4/10/2017 8:50 AM	DS_STORE File
ANSR 20180315.tar	3/16/2018 7:51 AM	TAR File
D005598v1.9 Release Notes, Complete So...	4/10/2017 8:50 AM	Adobe Acroba

Instrument Software Update

13. To update the instrument application, save the .bin file (from desktop folder) to a USB drive

    NOTE: *Do not rename .bin files. Ensure the .bin file is the only file on the USB drive.*

14. Insert the USB drive with the .bin file into the instrument and turn the instrument on. Follow the on-screen instructions on the instrument touch screen to update the instrument application.

15. The screen will display **Verifying update**, **and then you must touch the screen to update.**

    NOTE: *You only have 10 seconds to touch the screen). After installation is completed, screen prompts user to remove USB and touch the screen to continue. This completes the instrument application update.*

Instrument Settings

16. To set the date and time on the instrument, start at the instrument's home screen. Touch the **Settings** icon and press the down arrow to scroll down to **Administrator Settings**. The password is "admin." From here, the path is > ✓ (green check mark) > **Date and Time** > Set the correct date and time > ✓ (green check mark). The home icon will return the instrument to the home screen.

17. The audio and brightness on the instrument can be adjusted under the **Settings** icon. The latest software version can be found under **Settings** > About the device.

18. To control the reader from the PC software, set the instrument in **Remote Mode** (computer driven mode). From the home screen on the instrument touch screen, touch the **Settings** icon > **Remote Mode.** To ensure the instrument restarts every time in **Remote Mode**, touch the box next to **Boot into Remote Mode** so an X shows up in the box. Remain at this screen, and it is what will appear each time the instrument is booted up.

North America
Neogen Headquarters
800/234-5333 (USA/Canada)
foodsafety@neogen.com
foodsafety.neogen.com/en

Europe, Middle East and Africa
Neogen Europe
+ 44 (0) 1292 525 600
info_uk@neogeneurope.com
foodsafety.neogen.com/uk

Mexico
Neogen Latinoamerica
+52 (55) 5254-8235
informacion@neogenlac.com
foodsafety.neogen.com/sp

Brazil
Neogen do Brasil
+55 19 3935.3727
info@neogendobrasil.com.br
foodsafety.neogen.com/pt

China
Neogen Bio-Scientific Technology
+86 21 6271 7013
info@neogenchina.com.cn
www.neogenchina.com.cn

India
Neogen Food and Animal Security
+91 484 2306598, 2301582
info@neogenindia.com
www.neogenindia.com

© Neogen Corporation, 2019. Neogen, and ANSR are a registered trademarks of Neogen Corporation.