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

800/234-5333
517/372-9200

Surface Sampling


1. Turn the instrument on. Make sure it's at the main screen.
2. Select the correct site to be tested.
3. Sample your site by "drawing" a 4" X 4" square and then with a back-and-forth motion within the square. Consistency is critical.

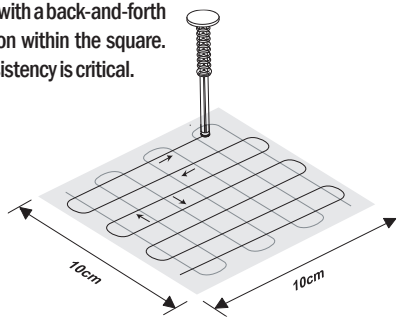
4. Hold the sampler vertically and fully depress to activate. Note: The sampler must be held vertically.

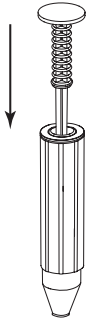
5. Mix for two seconds.


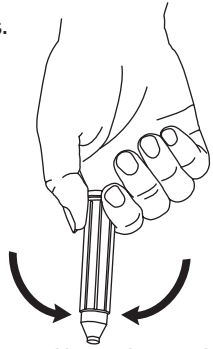
6. Press the eject button and insert the sampler into the instrument.

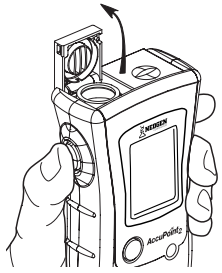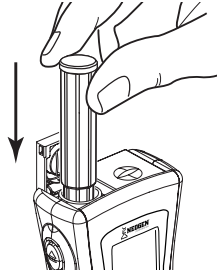
7. Depress the sampler all the way down and close the door.

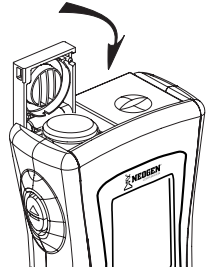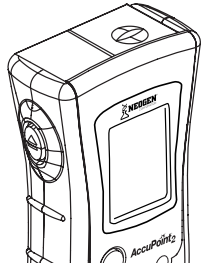8. The RLU score will display along with a symbol for pass, marginal or fail.

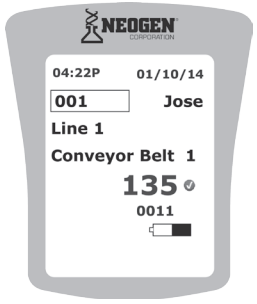AccuPoint Advanced

Sampling Quick Reference Card


1232A

NEOGEN CORPORATION


800/234-5333
517/372-9200

Access Sampler


1. Turn the instrument on. Make sure it's at the main screen.
2. Select the correct site to be tested.
3. Remove the Access sampler and take your sample.

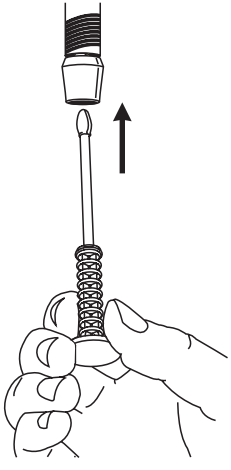
4. Hold the sampler vertically and fully depress to activate. Note: The sampler must be held vertically.

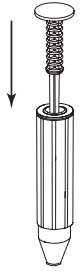5. Mix for two seconds.

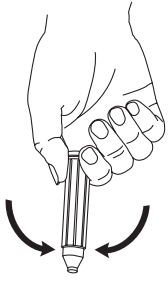6. Place in instrument immediately.

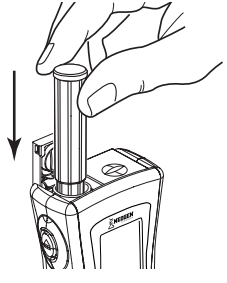Water Sampler

1. Turn the instrument on. Make sure it's at the main screen.
2. Select the correct site to be tested. A water droplet icon on the display indicates to use a water sampler.


3. Take your sample by dipping the water sampler into the liquid.

4. Hold the sampler vertically and fully depress to activate. Note: The sampler must be held vertically.

5. Mix for two seconds.

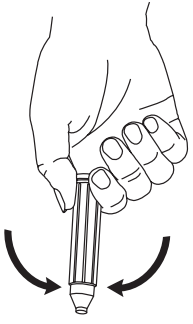6. Place in instrument immediately.