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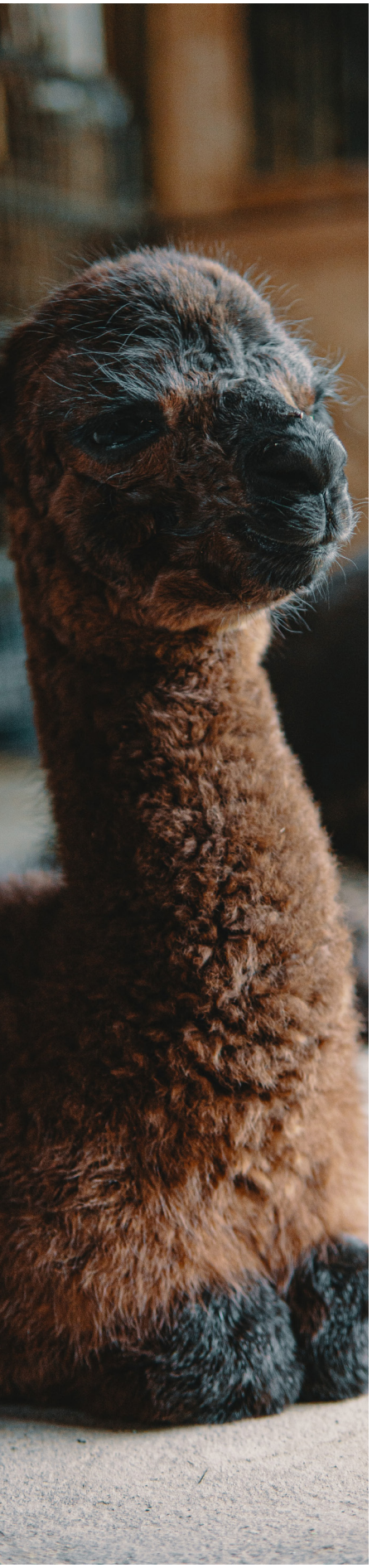
Alpaca Coat Colour Test

The Alpaca Coat Colour Test has been developed in collaboration with Dr. Kylie Munyard from Curtin University of Western Australia, who has conducted extensive research in the area of alpaca coat colour genetics.

Coat colour in Alpacas is a complex trait, involving two main genes responsible for base coat colour (AslP and MCIR), and an as yet unknown number involved with pattern. Alpaca fleece has 22 natural shades that ranges from black to white, grey, fawn to champagne. Breeding for a specific coat colour can be a complex process.

The 'classic grey' phenotype can be problematic in breeding due to its association with the blue eye white phenotype & associated possible health defects. Classic grey can be hidden or cryptic on white or light backgrounds.

With the release of the Alpaca Coat test, breeders have the opportunity to test their white or light fawn animals, those with uncertain patterns or mutations or animals they wish to determine the base coat colour to deduce common progeny colours. The test also identifies animals with 'cryptic grey' coat patterns that are generally too pale to see.

Example of Results

Breeders will receive a grey or non-grey status for tested animals, as well as a base coat phenotype for the following colours:

Base Colour Codes		
W	White	White fibre, Dark skin
PSW	Pink Skinned White	White fibre, Pink skin
F	Fawn	Fawn fibre, Dark skin
CF	Clear Fawn	Fawn fibre, Pink skin
BB	Bay/Brown	Red/Brown body fibre, Black fibre on extremities, Black skin
CH	Chestnut	Red/Brown fibre, Pink to Red/Brown
B	Black	Black fibre and skin

Collection Methods

There are 2 main methods of collection

Blood Cards

A small amount of blood collected onto a blood card. If using blood cards, it is suggested that collecting blood from the inside of the lower lip yields a more successful sample, than from the ear.

Tissue Sample Units (TSU)

TSUs would require the collection of a small ear notch from the center of the animal's ear, that self-seals into a collection tube. TSUs tend to be a more efficient and effective method of collection.

Contact Information

For further information, please contact NEOGEN Canada at 1.855.324.9774 or email us at neogencanada@neogen.com

NEOGEN.com

© NEOGEN Corporation, 2021. NEOGEN is a registered trademarks of NEOGEN Corporation, Lansing MI, 48912 U.S.